

## Product Description of SMMS Hydrophobic Waterproof Nonwoven For Medical Gown

Ours is SSMMS which is better quality than normal SMS. SSMMS would be better evenness, better tensile strength, better isolation ability than SMS.

but actually, for you, and most of the users, can not distinguish SMS and SSMMS by your eyes. because they look the same.

Normally, our 45gsm SSMMS are used the most for making surgical gown and drapes.

Our 45gsm can pass US standard AAMI level 3, and European standard EN 13795-1:2019 Surgical clothing and drapes - Requirements and test methods - Part 1: Surgical drapes and gowns the reports of AAMI and EN 13795 would be available.

now we are offering this SSMMS fabric to Turkey, USA, Italy, Peru, Azerbaijan, Columbia, Russia, United Kingdom.

**Product :** SMMS Hydrophobic Waterproof Nonwoven For Medical Gown

**Model No. :** GT-hysms45b

**Composition :** pp

**Technics :** spunbond + spunbond + meltblown + meltblown + spunbond

**Average weight of cloth :** 45 GSM

**Size :** 160 cm

**Length :** 1000 m

**Diameter of Each Roll :** 55 cm

**Weight :** approx 73 kg

**Volume :** approx 0.48 m<sup>3</sup>

**Colors :** Blue

**Standard :** US standard AAMI level 3, and European standard EN 13795-1:2019 Surgical clothing and drapes - Requirements and test methods - Part 1: Surgical drapes and gowns

**MOQ :** 1000 KG (customized)

**Package :** Packed in rolls, with 2 " or 3 " ID core inside, poly bags outside

**Sample :** Free sample with freight collect

**Characteristic :** anti-bacteria, anti-pull, anti-static, tear-resistant, breathable, eco-friendly, waterproof

**Applications :** surgical gowns, theatre scrubs, surgical drapes

our SSMMS products can reach the following standards of the surgical gown

我们的SSMMS产品能达到的以下手术衣标准

|   | Standard Test Items<br>标准测试项目      | AAMI PB70<br>美标 | EU 13795<br>欧标 | YY/T 0506<br>国标 |
|---|------------------------------------|-----------------|----------------|-----------------|
| Protection Performance<br>防护性能                | Hydrostatic Pressure Test<br>静水压测试 | √               | √              | √               |
|   | Liquid Permeability<br>液体渗透性       | √               | √              | √               |
|   | Blood Permeability<br>血液渗透性        | √#              | ○              | ○               |
|   | Microbial Permeability<br>微生物渗透性   | √#              | √              | √               |
| Material Strength<br>材料强度                     | Tensile Strength<br>拉力强度           | ○               | √              | √               |
|   | Blasting Strength<br>爆破强度          | ○               | √              | √               |
| The Comfort Level<br>舒适度                      | Breathability<br>透气性               | ○               | ○              | ○               |
|   | Moisture Permeability<br>透湿性       | ○               |                | ○               |
| Microbial Transmission Performance<br>微生物传递性能 | Deflocculation<br>脱絮               | ○               | √              | √               |
|   | Particle Cleanliness<br>颗粒洁净度      | ○               | √              | √               |
|   | Microbial Cleanliness<br>微生物洁净度    | ○               | √              | √               |

√: Test items included in the standard;

√: 标准中包含该测试项目;

#: Only applicable to Level 4 protection products;

#: 仅适用于四级防护产品;

○: Test items are not included in the standard

○: 标准中不包含改测试项目

美国医疗器械促进协会 (AAMI) 对于医用纺织产品如：手术衣、洞巾、隔离衣等液态阻隔性能等级分为四级  
 American Association for the Advancement of Medical Devices (AAMI) divides the liquid barrier performance into four levels for the medical textile products, such as surgical gown, drape and isolation clothing, etc.



第四级 (Level 4)  
 必须通过血液与病毒渗透两项试验, 抗静水压试验必须大于140cmH<sub>2</sub>O;  
 Must pass the test of Blood and virus penetration, hydrostatic pressure test must be more than 140cmH<sub>2</sub>O;

| 测试标准 Test         | 测试结果 Results | 测试要求 AQL Requirements |
|-------------------|--------------|-----------------------|
| ASTM F1671: 2003* | PASS         | 4%                    |
| ASTM F1670: 2003* | PASS         | 4%                    |

高防护



第三级 (Level 3)  
 疏水性渗透量必须小于1.0g, 抗静水压试验必须大于50cmH<sub>2</sub>O;  
 The hydrophobic permeability must be less than 1.0g, hydrostatic pressure test must be more than 50 cmH<sub>2</sub>O

| 测试标准 Test      | 测试结果 Results | 测试要求 AQL Requirements |
|----------------|--------------|-----------------------|
| AATCC 42:2000  | <1.0g        | 4%                    |
| AATCC 127:1998 | >50cm        | 4%                    |

加强防护



第二级 (Level 2)  
 疏水性渗透量必须小于1.0g, 抗静水压试验必须大于20cmH<sub>2</sub>O;  
 The hydrophobic permeability must be less than 1.0g, hydrostatic pressure test must be more than 20cmH<sub>2</sub>O

| 测试标准 Test      | 测试结果 Results | 测试要求 AQL Requirements |
|----------------|--------------|-----------------------|
| AATCC 42:2000  | <1.0g        | 4%                    |
| AATCC 127:1998 | >20cm        | 4%                    |

标准防护




第一级 (Level 1)  
 医用纺织品必须经过疏水性 (抗冲击渗透) 试验, 渗透量必须小于4.5g  
 Medical textiles must pass the hydrophobicity (impact penetration resistance) test, and the penetration amount must be less than 4.5g

| 测试标准 Test     | 测试结果 Results | 测试要求 AQL Requirements |
|---------------|--------------|-----------------------|
| AATCC 42:2000 | <4.5g        | 4%                    |

标准防护



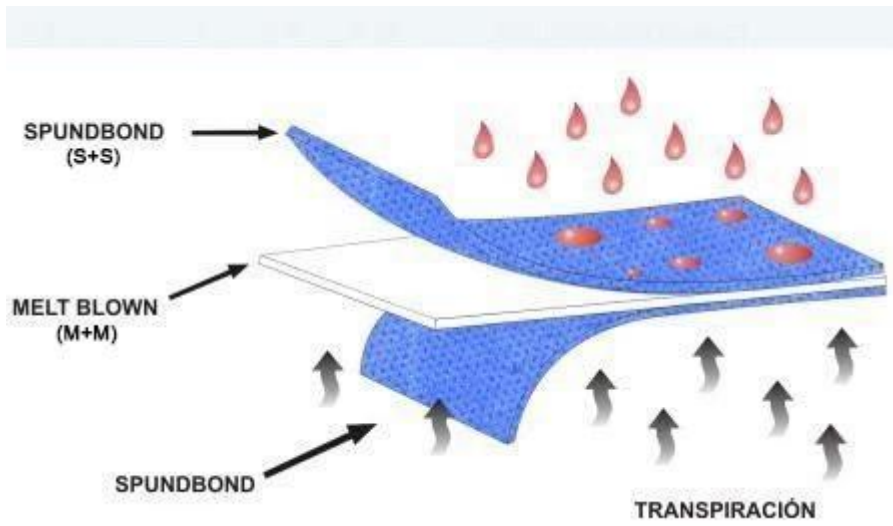


Product Show of SMMS Hydrophobic Waterproof Nonwoven For Medical Gown

[Hydrophobic SMS Supplier](#)





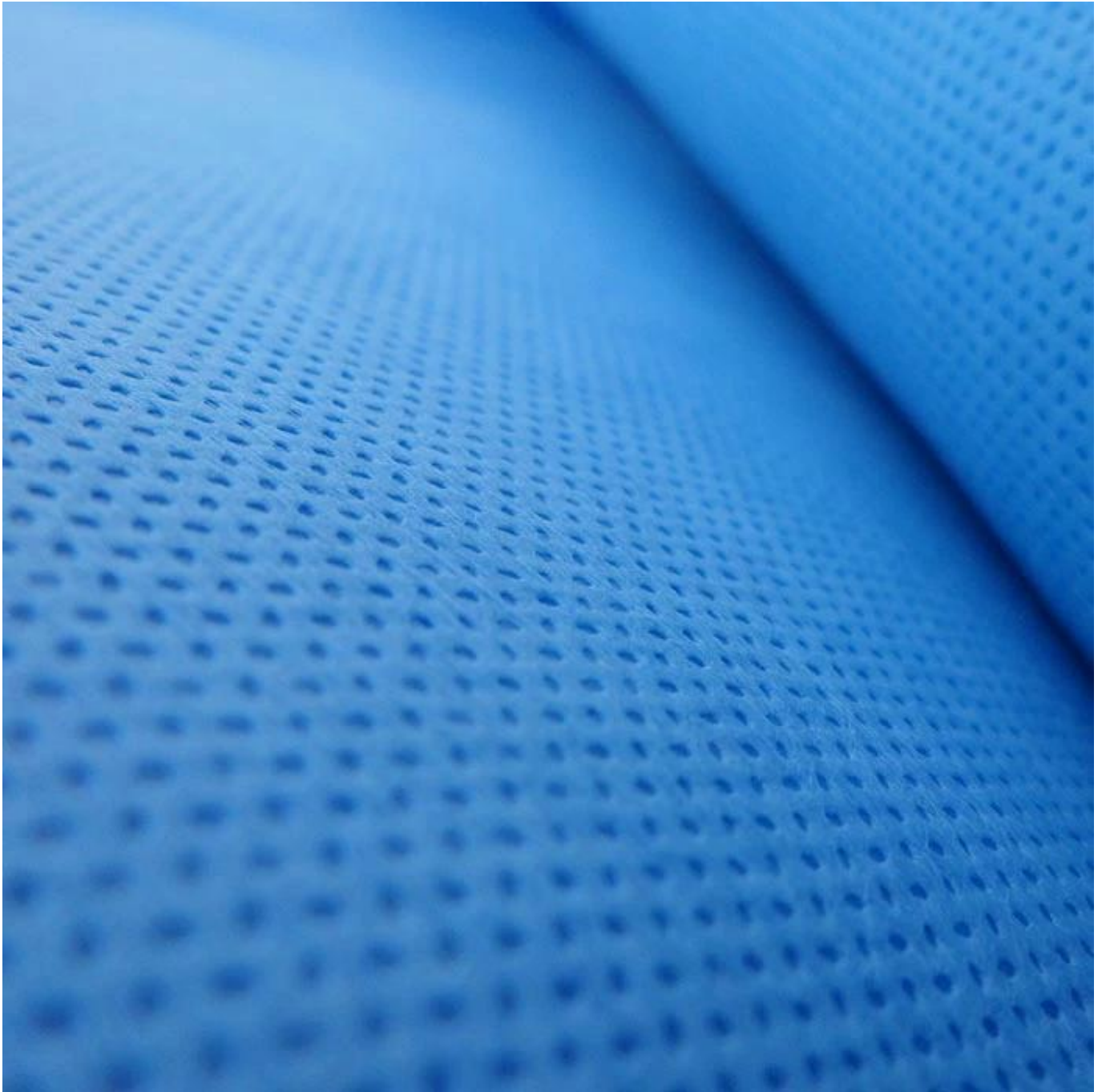

Applications of SMMS Hydrophobic Waterproof Nonwoven For Medical Gown

[Nonwoven Sms Manufacturer](#)



**Disposable Underwear**



**Nonwoven Shoes Cover**



**Protective Clothing**



**Nonwoven Pillow bag**



**Nonwoven Bed Sheet**



**Face Mask**



**Medical Wrapping**



**Disposable Bouffant Cap**



**Nonwoven Sleeve**

Product Packaging and Shipping of SMMS Hydrophobic Waterproof Nonwoven For Medical Gown





Production Line



Product Packaging



Loading

Certificate of SSMS Nonwoven Fabric



|                              |  |  |
|------------------------------|--|--|
| <b>检测报告(Testing Report)</b>  | 防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a><br>防伪码(Security code): 0561010345 |  |
| 报告编号(Report No): ZFLJ2626978 | Page 1 of 4  |  |

**Applicant Information (客户信息)**

Applicant Name(委托单位) : FOSHAN GUIDE TEXTILE CO., LTD(佛山凯缙纺织品有限公司)  
Applicant Address(委托单位地址) : No.8, Hezhong RD, Yanbu, Dali Town, Nanhai Dist., Foshan, G.D.528247  
CHINA (佛山市南海区大沥镇盐步河中路8号)

**Sample Information (样品信息)**

Sample Description(样品描述) : SMS non-woven fabric 45g (SMS 无纺布 45 克)  
Style No.(款号) : KHSSMMS-45G  
Sample Quantity(样品数量) : 8 m (8 米)  
Color and Description(颜色及描述) : Medical blue: 45g (医用蓝色: 45 克)  
Size specification(号型规格) : KHSSMMS-45G

- Sample Receiving date (收到样品的日期) : 2020-06-02  
- Report date (报告日期) : 2020-06-05  
- The original sample is stucked on the last paper. (送检样品原样粘贴在本报告的末页)

**Test Performed (检验标准)**

Judgement according to(评定依据):  
GB 19082-2009 Technical requirements for single-use protective clothing for medical use  
医用一次性防护服技术要求  
GB/T 38462-2020 Textiles-Nonwoven fabrics for isolation gowns Level II  
纺织品 隔离衣用非织造布 II 级  
- Selected test(s) as requested by applicants. (本实验室根据客户要求完成以下检测内容)

**Pronounce (声明)**

The results shown in this report refer only to sample(s) tested unless otherwise stated.All the items are performed in the standard conditions, except the noted cases.(除非特别说明, 测试结果只对来样负责。所有测试项目均在标准规定的环境下进行, 有注明除外。)  
Except for the requirement of the client, the test results and the conformity judgement of this report do not take the uncertainty of the test results into account.(除非客户要求, 本报告检测结果及符合性判定不考虑测量结果的不确定度。)

Signed for and on behalf of  
CNTAC Testing Service Co.,Ltd.(Foshan)

Approved by  
(批准)

张志强





|                              |   |  |
|------------------------------|---|--|
| <b>检测报告(Testing Report)</b>  | 防伪查询网址(Security website): www.fcl-sz.org.cn<br>防伪码(Security code): 0561010345 |  |
| 报告编号(Report No): ZFLJ2626978 | Page 2 of 4   |  |

-- Test Result (测试结果) --

| Test Result(测试结果) | Requirements(评定条件) | Judgement(判定) |
|-------------------|--------------------|---------------|
| ✓ GB 19082-2009   |                    |               |

**1. Tensile Strength(断裂强力) GB/T 3923.1-1997 Strip test(条样法)**

5cm strip test, CRE tensile test machine, testing speed 100mm/min, gauge length 100mm.(5 厘米条样法, CRE 拉伸测试仪, 测试速度 100mm/min, 初始隔距 100mm)

Unit(单位): <N>

|                |     |      |          |
|----------------|-----|------|----------|
| Lengthwise(纵向) | 104 | ≥ 45 | Pass(合格) |
| Crosswise(横向)  | 51  |      |          |

**2. Breaking Elongation(断裂伸长率) GB/T 3923.1-1997 Strip test(条样法)**

5cm strip test, CRE tensile test machine, testing speed 100mm/min, gauge length 100mm.(5 厘米条样法, CRE 拉伸测试仪, 测试速度 100mm/min, 初始隔距 100mm)

Unit(单位): <%>

|                |    |      |          |
|----------------|----|------|----------|
| Lengthwise(纵向) | 94 | ≥ 15 | Pass(合格) |
| Crosswise(横向)  | 78 |      |          |

**3. Filtration Efficiency(过滤效率) GB 19082-2009 Section 5.7(条款 5.7)**

Air flow: 15L/min, Aerosol: NaCl(气体流量: 15L/min, 气溶胶颗粒: NaCl)

Unit(单位): <%>

|                    |      |      |          |
|--------------------|------|------|----------|
| Minimum value(最小值) | 72.6 | ≥ 70 | Pass(合格) |
|--------------------|------|------|----------|

**4. Hydrostatic Pressure(抗渗水性) GB/T 4744-1997**

Distilled water temperature is 20°C. The rate of increase of water pressure is 6.0kPa/min.

Testing surface: face(测试水温度 20°C, 水压上升速率 6.0kPa/min, 测试面: 正面)

Unit(单位): <kPa>

|                 |      |        |          |
|-----------------|------|--------|----------|
| Mean value(平均值) | 6.93 | ≥ 1.67 | Pass(合格) |
|-----------------|------|--------|----------|

**5. Resistance to Penetration by Synthetic Blood(抗合成血液穿透性) GB 19082-2009 Appendix A(附录 A)**

Synthetic blood surface tension: 0.040N/m(合成血液表面张力: 0.040N/m)

Unit(单位): <grade(级)>

The sample cannot be pressurized due to air leakage, so it is not suitable for testing

(样品漏气无法加压, 不适合测试)





|                              |  |               |
|------------------------------|--|---------------|
| <b>检测报告(Testing Report)</b>  | 防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a><br>防伪码(Security code): 0561010345 |               |
| 报告编号(Report No): ZFLJ2626978 | Page 3 of 4  |               |
| Test Result(测试结果)            | Requirements(评定条件)   | Judgement(判定) |

**6. Water Vapor Permeability(透湿量) GB/T 12704-1991 Method A(方法 A)**

Desiccant method at 38°C, Relative humidity at 90%, wind speed at 0.3m/s, Testing surface: back.(吸湿法, 测试仓温度: 38°C, 湿度: 90%, 风速: 0.3m/s, 测试面: 反面)

Unit(单位):  $<g/(m^2 \cdot d)>$

Result(测试结果)  $7.09 \times 10^3$   $\geq 2500$  Pass(合格)

**✓ GB/T 38462-2020**

**1. Hydrostatic Pressure(静水压) GB/T 24218.16-2017**

Distilled water temperature is 21°C. The rate of increase of water pressure is 1.0kPa/min.

Testing surface: face (测试水温度 21°C, 水压上升速率 1.0kPa/min, 测试面: 正面)

Unit(单位):  $<kPa>$


Minimum value (最小值) 6.5  $\geq 1.8$  Pass(合格)

Remark (备注): The hydrostatic pressure item is unauthorized for CNAS.  
(CNAS 未授权静水压项目。)





佛山中纺联检验技术服务有限公司  
 CNTAC Testing Service Co.,Ltd.(Foshan)

|                              |  |   |
|------------------------------|--|---|
| <b>检测报告(Testing Report)</b>  | 防伪查询网址(Security website): <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a><br>防伪码(Security code): 0561010345 |  |
| 报告编号(Report No): ZFLJ2626978 | Page 4 of 4  |   |

Original (送检样品原样)



=====  
 ===== End of Report (检测报告结束) =====  
 =====



广东省佛山市南海区西樵镇南方技术创新中心一、三、四楼 1, 3&4F., Nanfang Technology Innovation Center, Xiqiao, Nanhai District, Foshan, Guangdong, China  
 Tel: (86 757)86850633/86806656 Fax: 8685 0633 Http://www.fcl.org.cn

若对检测报告有异议, 应于收到报告之日起 15 日内向检测单位提出, 逾期不予受理。

## Our Company

**FOSHAN GUIDE TEXTILE CO., LTD.**, which is located in Foshan city. We are engaging in manufacture and sales of non woven fabric and non woven related products, since 2005. Up to now, we have

professional elite sales team, several advanced non woven production lines with width of 1.6M, 1.9M, 2.1M, 2.4M, 2.5M, 3.2M, 3.7M, 4.8M, etc. Our products are well sold both in domestic and overseas markets, widely used in areas of agriculture and gardening, medical and hygiene, personal care and beauty, industry application, packaging and wrapping, foundation and building, hotel and restaurant, labor protection and household, decoration, etc. Hot-sale products are pp (polypropylene) non woven fabric, pet (polyester) non woven fabric, sms non woven fabric, elastic non woven fabric, flower wrapping material, disposable underwear, disposable bed sheet, weed control fabric, disposable table cloth, fruit protection bag.



## FAQ of SMMS Hydrophobic Waterproof Nonwoven For Medical Gown

| Question   | Answer  |
|--|---|
| Q1: Do you have it in stock ?  | A1: Yes,we have.  |
| Q2: What is your MOQ ?   | A2: Customized:1ton   |
| Q3: How can I get the sample ?   | A3: If you need sample to test, we can make it as per your request. But we need to charge sample fees. If it is our regular product in stock, you just pay freight cost and sample is free.                       |
| Q4: What is your package detail ?  | A4: individually packed for each pcs,or customer's requirements   |
| Q5: What should I do if I cannot find what I looking for on your website ? | A5: Please feel free to contact us by email and send detailed information for the product you look for.We will check whether it is available.   |
| Q6: What about lead time for mass production ?                             | A6: Honestly, it depends on the order quantity and the season you place order. Generally speaking,the lead time is about 2-4 weeks. So we suggest that you start inquiry as earlier as possible.                  |
| Q7: Is it all right to make customer's own brand name ?                    | A7: Sure. We warmly welcome the OEM clients. We can discuss the package, design, etc.   |
| Q8: When can I get the price ?   | A8: Usually we quote within 24 hours after we get your inquiry detail. If you are very urgent to get the price, please call us or tell us in your email so that we will quote your the price as soon as possible. |